


## Oxygen Monitor Module Datasheet (Part Number OGI-A-1008)

The Oxygen Monitor module is an optional accessory for the OGI Microbioreactor Mk2 Part number OGI-A-1006. This module requires culture flasks with oxygen sensor spots (Part number A200017). After calibration with reference samples, the module can be used to measure the amount of dissolved oxygen in the culture media at regular intervals during an experiment.

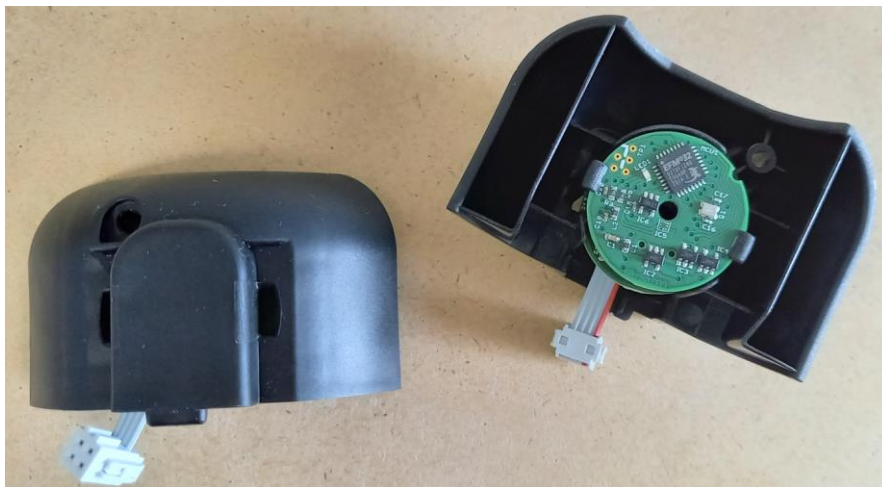

### Technical Specification

- Power Requirements:
  - o This unit is powered from OGI Microbioreactor Mk2 Part number OGI-A-1006
- The device is intended for indoor use only
- Optimal Performance at:
  - o Operating temperature: 15 to 40°C
  - o Relative humidity: 30% - 70%
  - o Operating altitude up to 2000 meters